

How to Insert a Respimat[®] Cartridge



Interior Health
Every person matters

Patient Information


www.interiorhealth.ca



1. Keep the cap closed. Push down on the safety catch and pull the clear base off the Respimat[®].



2. Take the cartridge out of the box. Slide it into the top part of the Respimat[®].

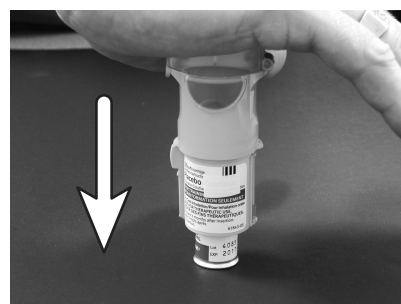

3. Put the Respimat[®] and cartridge on a table or counter. Press down firmly on the top of the Respimat[®] so the cartridge goes all the way in.

Note: you will see some of the cartridge at the bottom of the device. This is correct.




4. Slide the clear base back on the Respimat[®]. It will click into place.

Do not remove the base or the cartridge once you have it secure.



5. With the cap closed; turn the base of the Respimat[®], in the direction of the arrows, a half turn. You will hear it click.

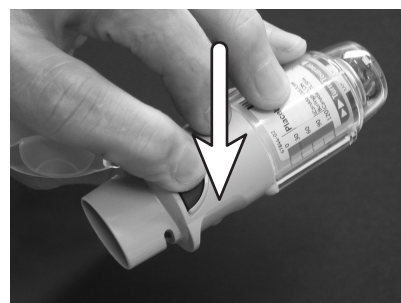

6. Open the cap and point the Respimat[®] down. Press the button. Repeat (steps 5 to 6) three more times. The Respimat[®] is ready to use.

This document is for informational purposes only. For detailed instructions on use of the Respimat[®] and information about your medication please see the package insert or ask your pharmacist.

Healthlink BC Call 8-1-1 www.healthlinkbc.ca	Nurse	24 hours a day	Daily	Call 8-1-1 to speak with a nurse, ask a dietitian about nutrition, or a pharmacist about your medication.
	Dietitian	9 am–5 pm	Mon–Fri	
	Pharmacist	5 pm–9 am	Daily	
	Hearing Impaired	Call 7-1-1		